

pTT5

PVT1032 2ug

Lot no 1901093

Store at -20°C

Expiry date

pTT5 Information

Promoter: CMV

Replicon: PMB1 ori

Terminator: beta globin poly (A) signal

Plasmid classification: mammalian cell carrier; protein over expression plasmid

Plasmid size: 4401bp

Prokaryotic resistance: ampicillin Amp

Clone strain: DH5 alpha

Culture conditions: 37, aerobic, LB


Expression host: mammalian cells

Culture conditions: 37, aerobic, LB

Primers for 5'sequencing: CMV-F:CGCAAATGGGCGGTAGGCGTG

Primers for 3'sequencing: Based on sequence design

pTT5 sequence Map

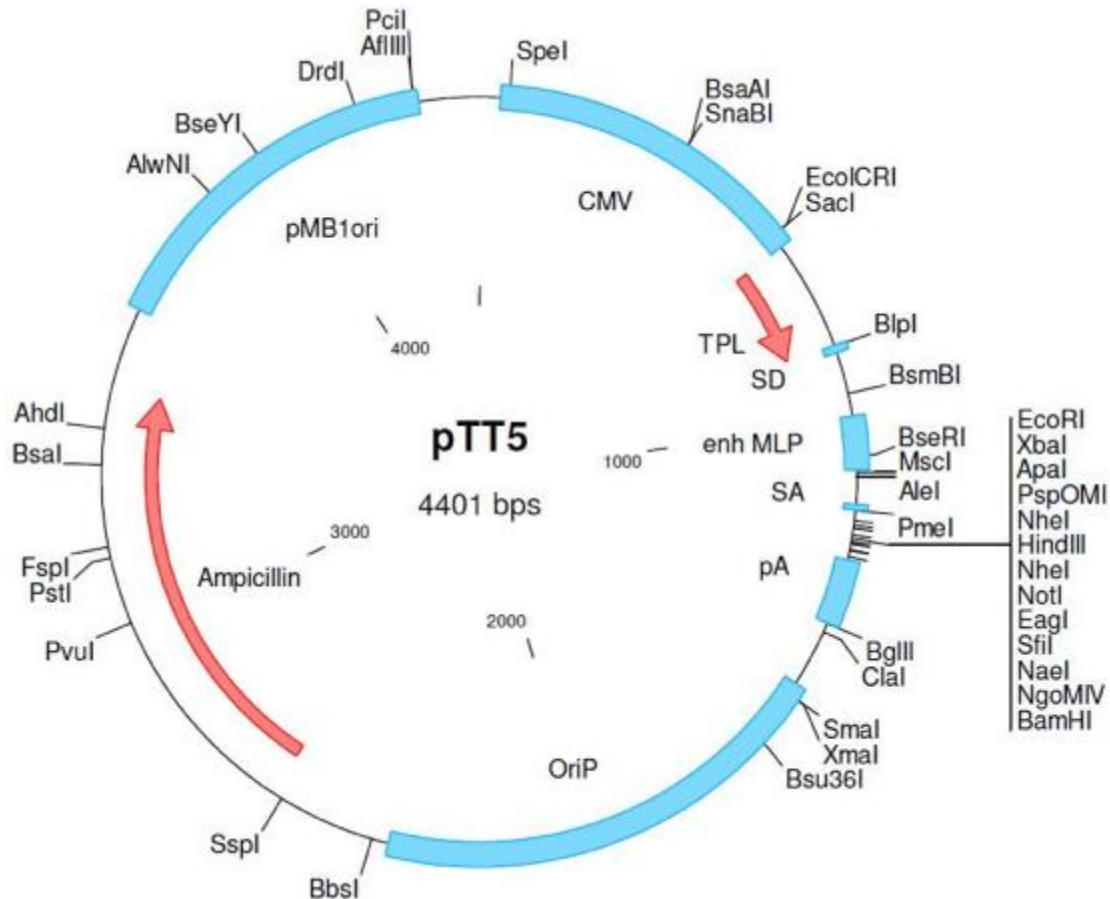

pTT5 sequence

LOCUS Exported 4401 bp ds-DNA circular SYN 08-July-2015

DEFINITION synthetic circular DNA

ACCESSION .

VERSION .

KEYWORDS pTT5

SOURCE synthetic DNA construct

ORGANISM synthetic DNA construct

REFERENCE 1 (bases 1 to 4401)

AUTHORS Huanghua

TITLE Direct Submission

JOURNAL Exported 2016-9-18 from SnapGene Viewer 2.8.1

FEATURES Location/Qualifiers

source 1..4401

/organism="synthetic DNA construct"

/mol_type="other DNA"

enhancer 42..421

/note="CMV enhancer"

Nova lifetech Limited.

Add: MNJ 2059 RM 1062 16/F

OMEGA PLAZA 32 DUNDAS ST MONGKOK KI

HONGKONG

Tel:(852)81700345 Email: sales@lifescience-market.com

```

        /note="human cytomegalovirus immediate early enhancer"
promoter    422..625
        /note="CMV promoter"
        /note="human cytomegalovirus (CMV) immediate early
        promoter"
polyA_site   1277..1332
        /note="beta-globin poly(A) signal"
        /note="rabbit beta-globin polyadenylation signal"
promoter    2501..2605
        /gene="bla"
        /note="AmpR promoter"
CDS         2606..3466
        /codon_start=1
        /gene="bla"
        /product="beta-lactamase"
        /note="AmpR"
        /note="confers resistance to ampicillin, carbenicillin, and
        related antibiotics"
        /translation="MSIQHFRVALIPFFAAFCLPVFAHPETLVKVKDAEDQLGARVGYI
        ELDLNSGKILESFRPEERFPMMSSTFKVLLCGAVLSRVDAGQEQLGRRIHYSQNDL
VEYS
        PVTEKHLTDGMTVRELCSAAITMSDNTAANLLLTIGGPKELTAFLHNMGDHVT
RLDRW
        EPELNEAIPNDERDTTTPAAMATTLRKLTTGELLTLASRQQLIDWMEADKVAGP
LLRSA
        LPAGWFIADKSGAGERGSRGIIAALGPDGKPSRIVVIYTTGSQATMDERNRQIAEI
GAS
        LIKHW"
rep_origin   3637..4225
        /direction=RIGHT
        /note="ori"
        /note="high-copy-number ColE1/pMB1/pBR322/pUC origin of
        replication"
ORIGIN
    1 gtacattat attggctcat gtccaatag accgcatgt tgacattgat tattgactag
    61 ttattaatag taatcaatta cggggtcatt agttcatagc ccatatatgg agttccgcgt
    121 tacataactt acggtaaag gcccgcctgg ctgaccgcc aacgacccc gccattgac
    181 gtcaataatg acgtatgtc ccatagtaac gccaataggg actttccatt gacgtcaatg
    241 ggtggagtat ttacggtaa ctgccactt ggcagtacat caagtgtatc atatgccaag
    301 tccgcccct attgacgtca atgacggtaa atggccgcc tggcattatg cccagtacat
    361 gacctacgg gactttccta cttggcagta catctacgta ttagtcatcg ctattaccat
    421 ggtgatgagg tttggcagt acaccaatgg gcgtggatag cggtttgact cacggggatt

```

481 tccaagtctc cacccattg acgtcaatgg gagtttgttt tggcaccaa atcaacggga
541 ctttcaaaaa tgtcgtata accccgcccc gttgacgcaa atgggcgga ggcgtgtacg
601 gtgggagtc tatataagca gagctcgttt agtgaaccgt cagatcctca ctctctccg
661 catcgctgtc tgcgagggcc agctgttggg ctgcggttg aggacaaact ctctcggtc
721 tttccagtac tcttggatcg gaaaccgctc ggctccgaa cggactccg ccaccgaggg
781 acctgagcga gtccgcatcg accggatcgg aaaacctctc gagaaaggcg tctaaccagt
841 cacagtgcga aggtaggctg agcaccgtgg cgggcggcag cgggtggcgg tcggggttgt
901 tttggcgga ggtgctgctg atgatgtaat taaagtaggc ggtcttgaga cggcggatgg
961 tcgaggtgag gtgtggcagg cttgagatcc agctgttggg gtgagtact ccttcaaaa
1021 gcgggcatta ctctcgcct aagattgtca gtttcaaaa acgaggagga tttgatattc
1081 acctggccc atctggccat aacttgagt gacaatgaca tccactttgc ctttctctc
1141 acaggtgtcc actcccaggt ccaagtttaa cggatctct agcgaattcc ctctagaggg
1201 cccgtttctg ctagcaagct tgtagcggc cgctcgaggc cggcaaggcc ggatccccg
1261 acctgacct ctggctaata aaggaaattt atttcattg caatagtgtg ttggaattt
1321 ttgtgtctc cactcgaag gacatatggg agggcaaac atttggtcga gatccctcg
1381 agatctctag cttaggagc gatccccgc cggacgaac taaacctgac tacgacatct
1441 ctgcccctc ttcgcgggc agtgcagtga atccctcag ttggttgga caacttgcca
1501 actgaacct aaacgggtag catatgctc cgggtagta gtatatacta tccagactaa
1561 ccctaattca atagcatatg ttaccaacg ggaagcatat gctatcgaat tagggtagt
1621 aaaaggtcc taaggaacag cgatgtaggg gggcgggcca agatagggg cggattgctg
1681 cgatctggag gacaattac acacactgc gcctgagcgc caagcacagg gttgttggc
1741 ctcatattca cgaggtcgt gagagcacgg tgggctaag ttgcatggg tagcatatac
1801 taccxaaata tctggtagc atatgctatc ctaatctata tctggtagc ataggctatc
1861 ctaatctata tctggtagc atatgctatc ctaatctata tctggtagt atatgctatc
1921 ctaattata tctggtagc ataggctatc ctaatctata tctggtagc atatgctatc
1981 ctaatctata tctggtagt atatgctatc ctaatctata tccggtagc atatgctatc
2041 ctaatagaga ttaggtagt atatgctatc ctaattata tctggtagc atatactacc
2101 caaatatctg gatagcatat gctatcctaa tctatatctg gtagcatat gctatcctaa
2161 tctatatctg gtagcatag gctatcctaa tctatatctg gtagcatat gctatcctaa
2221 tctatatctg gtagtatat gctatcctaa tttatatctg gtagcatag gctatcctaa
2281 tctatatctg gtagcatat gctatcctaa tctatatctg gtagtatat gctatcctaa
2341 tctgatccg gtagcatat gctatcctca tgataagctg tcaaacatga gaattaattc
2401 ttgaagacga aaggcctcg tgatacgct attttatag gtaagtca tgataaat
2461 ggtttcttag acgtcaggtg gcacttttc gggaaatgtg cgcggaacc ctatttgtt
2521 attttctaa atacattcaa atagtatcc gctatgaga caataacct gataaatgct
2581 tcaataat tgaaaagga agagtagg tattcaacat ttccgtgctg ccctattcc
2641 ctttttgcg gcattttgc ttctgttt tgctcacca gaaacgctgg tgaaagtaa
2701 agatgctgaa gatcagttg gtgcagagt gggttacatc gaactggatc tcaacagcg
2761 taagatcct gagagtttc gcccgaaga acgtttcca atgatgagca cttttaaagt
2821 tctgctatgt ggcggtat tatccgtgt tgacccggg caagagcaac tcggtcggc
2881 catacactat tctcagaatg acttggttga gtactacca gtcacagaaa agcatctac
2941 ggatggcatg acagtaagag aattatgacg tgctgccata accatgagtg ataactg

3001 ggccaactta cttctgacaa cgatcggagg accgaaggag ctaaccgctt tttgacaa
3061 catgggggat catgtaactc gccttgatcg ttgggaaccg gagctgaatg aagccatacc
3121 aaacgacgag cgtgacacca cgatgcctgc agcaatggca acaacgttc gcaactatt
3181 aactggcgaa ctacttactc tagcttcccg gcaacaatta atagactgga tggaggcgga
3241 taaagttgca ggaccacttc tgcgctcggc ccttccggct ggctggttta ttgctgataa
3301 atctggagcc ggtgagcgtg ggtctcgcgg tatcattgca gactggggc cagatggtaa
3361 gccctccgt atcgtagtta tctacacgac ggggagtcag gcaactatgg atgaacgaaa
3421 tagacagatc gctgagatag gtgcctcact gattaagcat tgtaactgt cagaccaagt
3481 ttactcatat atactttaga ttgattaaa acttcatttt taatttaaaa ggatctaggt
3541 gaagatcctt tttgataatc tcatgaccaa aatcccttaa cgtgagtttt cgttccactg
3601 agcgtcagac cccgtagaaa agatcaaagg atcttctga gatcctttt ttctgcgct
3661 aatctgctgc ttgcaaaaa aaaaccacc gctaccagcg gtggtttgtt tgccggatca
3721 agagctacca actcttttc cgaaggtaac tggcttcagc agagcgcaga taccaaatac
3781 tgttcttcta gtgtagccgt agttaggcca ccactcaag aactctgtag caccgcctac
3841 atacctcgt ctgtaatcc tgttaccagt ggctgctgcc agtggcgata agtcgtgtct
3901 taccgggttg gactcaagac gatagttacc ggataaggcg cagcggtcgg gctgaacggg
3961 gggttcgtgc acacagccca gcttgagcg aacgacctac accgaactga gatactaca
4021 gcgtgagcta tgagaaagcg ccacgcttc cgaagggaga aaggcggaca ggtatccggt
4081 aagcggcagg gtcggaacag gagagcgcac gagggagctt ccagggggaa acgcctggt
4141 tctttatagt cctgtcgggt ttcgccacct ctgacttgag cgctgatttt tgtgatgctc
4201 gtcagggggg cggagcctat ggaaaaacgc cagcaacgcg gccttttac ggttctggc
4261 cttttgctgg cttttgctc acatgttctt tctgctgta tcccctgatt ctgtggataa
4321 ccgtattacc gccttgagt gagctgatac cgctcggcgc agccgaacga ccgagcgcag
4381 cgagtcagtg agcgaggaag c

//

Caution:

Product is for research use only!